

Titanium Nitride Coated Implant Abutments: From Technical Aspects to Clinical Applications. A Literature Review

Rafael del Castillo, DDS, FACP ¹,¹ Konstantinos Chochlidakis, DDS, MS, FACP ¹,¹
Pablo Galindo-Moreno, DDS, PhD,² & Carlo Ercoli, CDT, DDS, FACP, MBA ^{1,3}

¹Department of Prosthodontics, Eastman Institute for Oral Health, University of Rochester, Rochester, NY

²Department of Oral Surgery and Implant Dentistry, School of Dentistry, University of Granada, Granada, Spain

³Department of Prosthodontics and Periodontology, Eastman Institute for Oral Health, University of Rochester, Rochester, NY

Keywords

CAD/CAM; dental implants; hard coatings; implant esthetics.

Correspondence

Rafael del Castillo DDS, FACP, Department of Prosthodontics, Eastman Institute for Oral Health, University of Rochester, 625 Elmwood Ave, Rochester, NY 14623. E-mail: cashier@infomed.es

Disclosure statement: The authors declare no conflict of interest.

Accepted October 26, 2021

doi: 10.1111/jopr.13446

Abstract

Purpose: To review the most up to date scientific evidence concerning the technical implications and clinical applications derived from the use of titanium nitride hard thin film coatings on titanium alloy implant abutments.

Materials and Methods: A review was performed to answer the following focused question: “What is the clinical reliability of nitride coated titanium alloy abutments?”. A MEDLINE search between 1980 and 2021 was performed for investigations pertaining the clinical use of nitride coated titanium alloy implant abutments (TiN) in case reports, case series and short- and long-term non/randomized controlled clinical trials. Literature analysis led to additionally evaluate research related to the technical and biological aspects, as well as the physicochemical characteristics of TiN hard thin film coatings and their impact on titanium abutment biocompatibility, mechanical properties, macroscopic surface topography, and optical properties. Therefore, preclinical data from biomechanical and in vitro investigations were also considered as inclusion criteria.

Results: The limited number of clinical investigations published made a systematic review and meta-analysis not possible, therefore a narrative review was conducted. TiN coatings have been applied to dental materials and instruments to improve their clinical longevity. Implant abutments are coated with titanium nitride to mask the titanium oxide surface and enhance its surface characteristics providing the TiN abutment surface with a low friction coefficient and a very high chemical inertness. TiN coating is suggested to reduce early bacterial colonization and biofilm formation and enhance fibroblast cell proliferation, attachment and adhesion when compared to Ti controls. Additionally, studies indicate that hard thin film coatings enhance the mechanical properties (hardness and wear resistance) of titanium alloy and appears as a yellow color when deposited on the titanium alloy substrate. To date, clinical investigations show that nitride coated titanium abutments provide promising short-term clinical outcomes.

Conclusions: Published research on nitride-coated abutments is still limited, however, the available biomedical research, mechanical engineering tests, in vitro investigations, and short-term clinical trials have, to date, reported promising mechanical, biological, and esthetic outcomes.

The use of titanium abutments has limitations especially in the esthetic zone due to peri-implant mucosa discoloration when a thin phenotype is present.^{1,2} In these cases, monolithic zirconia abutments have been recommended.³ However, in vitro studies showed that the fracture resistance of zirconia abutments is half of that reported for titanium abutments.⁴ To overcome the mechanical limitations of one-piece zirconia abutments while maintaining its esthetic advantages, hybrid

metal-zirconia abutments have been recommended,^{5,6} although its fabrication requires additional laboratory steps. Therefore, research indicates that no abutment material is capable of providing the ideal esthetic and mechanical properties for all clinical scenarios.

Stock metal abutments have been used in implant dentistry for cement-retained prostheses. However, the relative inability to customize these abutment margins to the surrounding

peri-implant mucosa may lead to the placement of excessively deep margins, incomplete cement removal, cement retention, and biologic complications, such as periimplantitis.^{7,8} CAD/CAM technology allows the fabrication of titanium abutments with customized margin positions and emergence profiles.^{9,10} While titanium biocompatibility has been sufficiently demonstrated,¹¹ titanium CAD/CAM abutments can be coated with hard titanium nitride thin film coating (TiN) in order to mask the titanium oxide layer, and improve the abutment surface characteristics while maintaining titanium biocompatibility. Preliminary reports indicate that TiN coating enhance fibroblast cell proliferation, attachment and differentiation, as well as reduce bacterial adherence and biofilm formation. Therefore, these abutments are being considered as alternatives to traditional gold, titanium, high-strength polycrystalline ceramic abutments, and hybrid abutments, although scientific data is still limited. This review describes the most up to date evidence on TiN-coated custom abutments, focusing on its enhanced soft tissue biocompatibility, identifying their most relevant biomechanical and esthetic properties, and describing their most relevant clinical applications.

Materials and methods

A MEDLINE search was performed for English language publications between 1980 and 2021, pertaining the clinical use of nitride coated titanium alloy implant abutments in case reports, case series and short- and long-term non/randomized controlled clinical trials. The following keywords and terms were used: titanium, nitride, titanium nitride, implant abutments, coating, soft tissue, color, and biocompatibility. A preliminary screening analyzed all titles and abstracts obtained. Those which not provided sufficient and/or pertinent information related to the objective of the present investigation were discarded. On the other hand, after careful review of titles and abstracts of the referenced publications, those considered to “possibly” add valuable information to the area of interest of the present review, full-text was obtained, printed, and reviewed. This protocol enabled broadening the scope of the review to include preclinical data from biomechanical investigations and animal studies.

For a coherent revision of the present research, data was organized in six main topics: titanium nitride as a coating material for medical and dental applications; biocompatibility; surface hardness and wear resistance; optical properties; limitations; and clinical applications of titanium nitride coated abutments.

Results

Titanium nitride as a coating material for medical and dental applications

Hard material coatings such as TiN have been used on cutting tools as well as on machine parts and decorative elements. In medicine, TiN has been used to coat orthopedic implants (hip and knee replacement prostheses) in order to reduce attrition coefficient, wear and corrosion, and increase surface hardness.^{12–14}

Several coating technologies have been proposed for the deposition of thin ceramic coatings (TiN and ZrO₂) on titanium surfaces, such as physical vapor deposition (PVD), thermal application, chemical vapor deposition (CVD),¹⁵ and Nd:YAG laser gas.^{16,17} These techniques allow for the thin TiN coating to be formed by the reaction of pure titanium (as a substrate) and highly reactive monoatomic nitrogen gas in the vapor phase.^{18,19} In an extensive review on the physical properties of TiN coatings, Al-Jabbari *et al*¹⁵ reported the coating layers to have a mean thickness between 1 and 5 μm with ranges between 0.06 and 10 μm . Qualitatively, TiN coatings may be deposited as a monolithic (single) uniform layer or in a more complex multilayered structure.¹⁵

In dentistry, TiN coating was introduced in the 1980s.^{20–22} Since then, TiN coatings have been applied to dental materials and instruments used in implant dentistry, orthodontics, endodontics, periodontics, and prosthodontics. More recently, dental implant abutments have been coated with TiN to couple the biomechanical properties of titanium with the requirements for esthetics in anterior implant cases.^{23–25} In addition, TiN coatings changed the implant/abutment’s surface chemical composition and mechanical properties thus reducing the interaction with biologic fluids and increasing the hardness and wear resistance of the metal substrate.^{15,26–28} Reports have also shown TiN as having very high chemical inertness and low friction coefficient.^{24,29}

Biocompatibility

Biocompatibility and cellular response to different abutment materials and surface modifications varies.^{30,31} Titanium, a reactive metal, is always covered by a titanium dioxide layer of approximately 5 to 20 nm thickness.⁹ Hard thin film coatings (TiN and ZrN) can mask the titanium oxide surface and modify the surface characteristics while keeping the material’s biocompatibility.³² Indeed, several studies showed that fibroblast cell proliferation, attachment, and adhesion on TiN samples is similar^{33–37} when compared to machined Ti controls. Other reports³⁸ showed that TiN coating deposited by Physical Vapor Deposition revealed an increase in fibroblast growth on zirconium nitride and titanium nitride when compared to polished titanium.

In addition, TiN coating is suggested to reduce early bacterial colonization and biofilm formation. Gröbner-Schreiber *et al*, both in *in vitro*³⁹ and *in vivo*⁴⁰ studies, showed significantly fewer adherent bacterial colonies on TiN and ZrN coated surfaces. These results were confirmed by Scarano *et al*⁴¹ who showed a significant reduction of bacterial adhesion on TiN-coated titanium implants when compared to control, uncoated ones. Bacterial adhesion can be enhanced by surface roughness and generally, a maximum Ra value of 0.2 μm is suggested in order to minimize plaque formation and enhance a stable soft tissue seal around the abutment.^{42–45} In this sense, TiN coating has been shown to be generally below this threshold, with documented Ra values of $0.159 \pm 0.040 \mu\text{m}$ and $0.217 \pm 0.025 \mu\text{m}$ when measured with confocal microscopy and white light interferometry, respectively.⁴⁶

While the effect of TiN coating on clinical peri-implant bone loss has been shown to be similar to other implant



Figure 1 The gold color that characterizes TiN abutments is due to a natural phenomenon known as “light interference.” Titanium nitride color characteristics, wavelength and purity, are very close to those of standard pure gold.

surfaces,⁴⁷ TiN coatings have been suggested to reduce bacterial adherence and biofilm formation without compromising bone-marrow mesenchymal cell adhesion, proliferation, and differentiation.^{48–50}

Surface hardness and wear resistance

Several studies have discussed the surface characteristics of titanium implant abutments and suggested that instruments and materials used during oral prophylaxis procedures and maintenance programs alter titanium alloy surfaces.^{51,52} Enhanced wear resistance and abrasive hardness are described as important advantages of TiN/ZrN film coatings and, therefore, are suggested to improve the surface properties of pure titanium. In vitro investigations have demonstrated an improved abrasion resistance of TiN samples when compared to uncoated titanium controls.^{53–55}

Investigations have also evaluated TiN coating hardness and determined TiN hardness to be within a range of 15 to 26 GPa.^{56,57} In addition, authors have estimated that nitride coatings can increase hardness of pure titanium by approximately 10 times (TiN Vickers hardness of 1300).⁵³

Optical properties

TiN, when applied to a titanium surface, produces a natural phenomenon known as “light interference” and appears as a yellow color⁵⁸ (Fig. 1). An in vitro investigation has shown that TiN coatings can be obtained in several range of compositions and colors depending on the stoichiometry (ratio between the relative quantities of substances taking part in a reaction or forming a compound, in this case, titanium and nitrogen).²⁵ From the colorimetric point of view, TiN color coordinates were calculated in terms of dominant wavelength (which determines the color hue) and purity (which determines the saturation of color). Lousa *et al*²⁵ found that a dominant wavelength of 580 to 582 nm was obtained for all the TiN samples, very close to that of pure gold reference spectrum (579 nm wavelength). Additionally, a purity (color saturation) of 0.84 was obtained, which, according to the authors, approaches the color purity corresponding to the reflectance of standard pure gold (0.91).

The optical properties of newly developed esthetic coatings were also evaluated by Pecnik *et al*.⁵⁹ Three different coating systems (Ti-ZrO₂, Ti-Al-ZrO₂, and Ti-Ag-ZrO₂) were deposited on titanium substrates of different roughness (polished, machined, and sand-blasted) interference colors (pink, yellow and white) and mucosa thickness (1 mm, 2 mm, 3 mm). Pink and yellow colors coatings showed positive and partially significant overall color (chroma) change (ΔE) in the green-red (Δa) and blue-yellow (Δb) axis, and were considered by the authors advantageous to obtain an improved peri-implant soft tissue appearance. Similar results were reported by Ishikawa-Nagai *et al*,⁶⁰ where light pink and light orange showed improved esthetic outcomes.

Limitations

Some investigators have outlined concerns when TiN hard coatings are applied to metal alloys used in dentistry. Studies have shown that when TiN film coatings are deposited by PVD and CVD (chemical vapor deposition), possible detachment of TiN/ZrO₂ coatings from the metal substrate surface of transmucosal implant abutments can occur.⁶¹ Pecnik *et al*⁵⁹ indicated that ZrO₂ thin films detached easier from polished Ti samples, and suggest, in accordance with Yang *et al*⁶² that thin film coatings would adhere better on machined or sand-blasted rather than on polished titanium surfaces. Another suggested possible explanation for this delamination was provided by Tong *et al*⁶³ who reported that this phenomenon can occur due to inadequate adhesive strength between the TiN film when deposited on pure titanium. The low surface hardness of titanium (microhardness of approximately 200 HV) was unable to support the hard TiN film coating (microhardness of more than 400 HV).

Although possible detachment of TiN surface layer from the metal substrate is considered a disadvantage, the evidence is not conclusive due to the limited number of clinical investigations. To overcome this limitation, Mengel *et al*⁶¹ suggested employing a low contact pressure when performing oral hygiene procedures on TiN coated abutments. In addition, Tong *et al*⁶³ suggest using alternative coating procedures such the Duplex Treatment Technique, which includes the use of plasma nitriding followed by the application of TiN coating to improve the surface hardness of titanium and increase the adherence of the film substrate.

Other reported limitations include the development of a localized adverse allergic reactions associated with the use of TiN coated transmucosal abutments. Although TiN coating is described as biochemically stable,⁶⁴ a clinical report shows TiN as being the possible cause of an allergic contact stomatitis in a patient treated with TiN implant abutments used to support a fixed dental prosthesis.⁶⁵ The authors reported that after abutment insertion, the patient complained of pain and discomfort, and the mucosa adjacent to the TiN abutments was erythematous. However, 1 month after TiN abutment removal and insertion of titanium healing abutments, mucosa appeared healthy and all signs and symptoms had disappeared. Patch tests disclosed a positive reaction to TiN coated abutments. To the best of the authors knowledge, this is the only report describing this type of clinical event.

Clinical applications of titanium nitride coated abutments

The available literature evaluating clinical outcomes associated with the use of TiN abutments is scarce and, to date, there are no long-term randomized clinical trials published. However, the present investigation shows that promising short-term clinical data has been reported when related specifically to the mechanical behavior and technical complications associated with the clinical use of TiN abutments. The loading capabilities of titanium alloy (ductility and tolerance to compressive and tensile stresses)⁴ together with the TiN coating layer surface properties^{15,27,28} suggest that the biomechanical behavior of the TiN abutments provides adequate short-term clinical outcomes. In a 3-year prospective clinical investigation, Ferrari *et al*⁶⁶ indicated TiN abutments can be predictably used to restore implant supported restorations in all anterior and posterior locations. The authors found a statistically significant higher success rate (100%) at 3 years for CAD/CAM TiN abutments when compared to one piece CAD/CAM zirconia abutments (82.2%).

Adequate short-term clinical outcomes have also been reported for TiN abutments when related to soft-tissue biocompatibility. Investigations found no significant differences in biologic outcomes, radiographic indices and short-term survival rates⁶⁷ among conventional titanium, TiN-coated and zirconia abutments.⁶⁸ A retrospective study evaluating the effect of TiN/Zr single CAD/CAM abutments on the presence of interproximal papilla in the maxillary anterior region, showed that both types of abutments represent a predictable treatment alternative and appear to improve the esthetic results and the papilla presence in the anterior maxilla.⁶⁹ When compared to conventional custom titanium abutments in the anterior maxilla, the same group⁷⁰ showed that CAD-CAM TiN and zirconia abutments also enhanced the presence of interproximal papilla.

The influence of abutment material color and peri-implant soft-tissue thickness is reported to significantly influence the esthetic outcome of implant-supported restorations.^{71,72} Zirconia abutments show a peri-implant soft tissue color closest to the soft tissues around natural teeth, although this difference has not been considered significant when compared to gold alloy abutments.⁷³ In addition, yellow and pink abutment coating colors are shown to improve the optical appearance and the lightness of the overlying mucosa.⁵⁹ Promising optical behavior has been reported with CAD/CAM TiN coated abutments and the characteristic gold color of TiN abutments is suggested to be advantageous for abutment soft tissue esthetic integration,⁷⁴ especially in anterior esthetic areas with a thin gingival phenotype (≤ 2 mm).⁷⁵

To date, the vast majority of studies reporting on the use of TiN abutments employ CAD/CAM technology for abutment design and manufacturing.^{66,69,70} The principal observed advantages of CAD/CAM manufactured abutments include comparable clinical outcomes to stock abutments,¹⁰ customization, precision of fit, durability and an enhanced formation of peri-implant soft tissue topography and emergence profile.⁷⁶

Case reports have also been published to illustrate the use of TiN abutments to restore single implants in anterior and posterior locations.⁷⁷ According to Ramsey *et al*⁷⁸ the main reasons

for choosing TiN abutments, especially in the esthetic zone, include optimal biocompatibility, superior fracture resistance when compared to one-piece zirconia ceramic abutments, adequate wear and corrosion resistance, reduced bacterial adhesion and excellent esthetic integration with the adjacent soft-tissues.

Taking into considerations the reduced number of clinical investigations but also the advantages of CAD/CAM abutments, the fact that TiN abutments are manufactured using this technology would suggest that custom TiN abutments could constitute a predictable alternative to stock/custom, metal/zirconia abutments. TiN abutments could be a clinically relevant alternative in mechanically challenging but esthetically demanding situations such excessive restorative vertical height (RVH),⁷⁹ areas of high occlusal loads, maxillary anterior areas where nonaxial loads are often observed, implant-supported restorations fabricated on narrow diameter implants⁸⁰ and anterior esthetic areas with thin gingival phenotype, where the use of one-piece zirconia abutments may, otherwise, be contraindicated (Figs 2 to 5).

Discussion

Although the use of TiN hard thin film coating technology on metal alloy substrates has been applied in the medical field to coat orthopedic implants (hip and knee replacement prostheses) for decades,^{12,19,81,82} its implementation in implant dentistry is more recent. The number of clinical trials published in peer review journals utilizing TiN implants/abutments is still limited. In this review, these reduced number of clinical investigations supposed a limitation and the review protocol only partially complied with reporting guidelines for the development of protocols of systematic reviews and meta-analysis (PRISMA-P).⁸³ Therefore, a systematic review and meta-analysis could not be conducted, which can be considered a shortcoming of the current review.

The heterogeneity of the topics analyzed, the design of the preclinical and laboratory publications included in the current review and, the subjectivity inherent when authors selectively evaluate and present the available literature, may add additional limitations to this review in the sense that it may be biased by selective reporting. On the other hand, manual searching and careful reading of referenced preclinical and laboratory publications, facilitated broadening the scope of the initial search strategy. This made the authors include several investigations on such diverse areas as clinical medical orthopedics, material surface coating technology and behavioral sciences of biomedical/dental materials. This protocol, besides bringing out relevant information on the technical aspects of TiN coated implant abutments, to certain extent, may have “compensated” the weaknesses and limitations of the present investigation. Nevertheless, reporting the data as objectively as possible, was a primary goal of the authors.

To the best of the author’s knowledge, the first review on titanium nitride coatings in clinical dentistry was published in 1992.²² More recently, Al-Jabbari *et al*¹⁵ updated the available evidence on TiN and nitrogen ion implanted coated dental materials, however the implementation of this technology in clinical implant dentistry is relatively recent. Therefore, as a

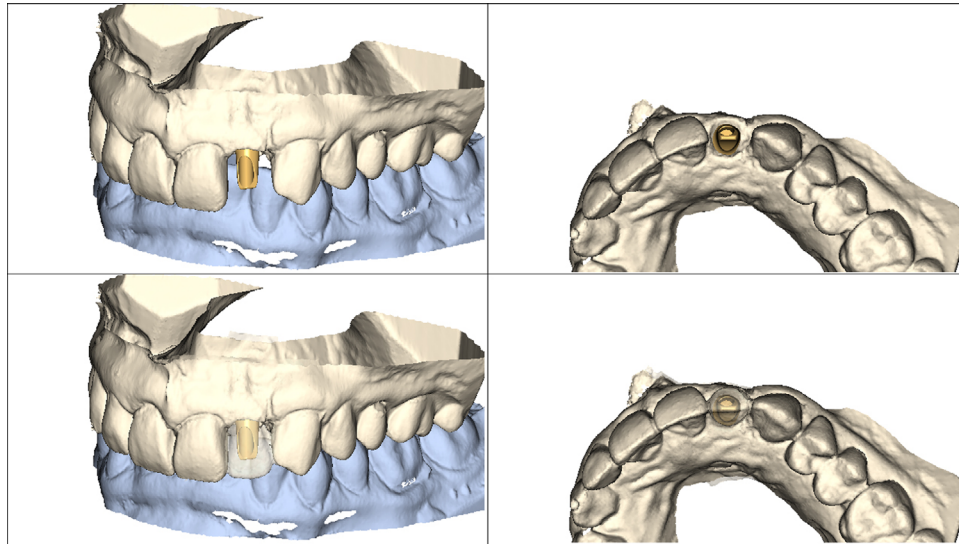


Figure 2 CAD-CAM technology is used to accurately design nitrided abutments based on patient's specific restorative requirements.

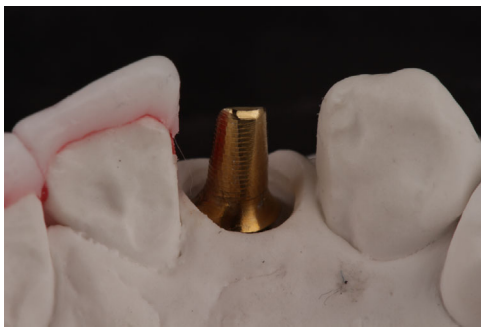


Figure 3 Abutment parameters such retention, resistance form, restorative margin location, core height and wall parallelism are optimized before definitive restoration is fabricated.



Figure 4 Clinical try-in of a custom TiN abutment (maxillary left lateral incisor) to ensure precision of fit prior to definitive restoration fabrication.

second objective, the present review is intended to complement the aforementioned publications by describing new evidence on the physicochemical characteristics, biomechanical properties and macroscopic surface topography of TiN hard thin film coatings and to highlight the potential applications of this



Figure 5 Adequate biocompatibility, fracture resistance and soft tissue esthetic integration suggests that nitrided abutments could represent a predictable alternative to more traditional abutment materials.

technology and their impact on titanium alloy implant abutment's biocompatibility, biomechanical properties and peri-implant soft tissue esthetic integration.

In light of the limited available evidence but also the potential applications of TiN hard thin film coating technology in clinical implant dentistry, there seems to be a need for additional research. At the preclinical level, research is necessary to elucidate the most ideal metal substrate (machined, sandblasted or polished titanium alloy) for hard thin film coating, the needed adhesive strength to prevent titanium nitride film detachment after long-term clinical use, and the clinical outcomes derived from these material combinations.

It is worth noting that the majority of recent clinical investigations using TiN coated abutments employ CAD/CAM technology for abutment design and manufacturing.^{66-70,74,77,78} Although the level of evidence is low due to the fact that few randomized controlled clinical trials are published, systematic reviews indicate that CAD/CAM abutments demonstrate overall good survival and success rates, and similar, if not superior, clinical outcomes when compared to conventional abutments.¹⁰ Therefore, at clinical level, additional

randomized controlled trials with long-term observational periods are necessary to confirm, not only, the promising short-term clinical outcomes of CAD/CAM nitride coated implant abutments, but also the advantages of custom abutments over stock abutments.

Conclusions

The combination of adequate clinical outcomes and reliability of titanium alloy abutments has led titanium to be the preferred material for custom implant abutment fabrication. Although published research and reported clinical outcomes on TiN prosthetic abutments is still limited, the available biomedical research as well as mechanical engineering tests, *in vitro* investigations, and short-term clinical trials have, to date, reported promising mechanical, biological and esthetic outcomes for TiN abutments. However, further *in vitro* investigations and long-term randomized clinical trials are needed to validate these initial observations, especially those investigating the most ideal metal substrate (machined, sand-blasted, or polished titanium alloy) for hard thin film coating, the needed adhesive strength to prevent titanium nitride film detachment, and the clinical outcomes derived from these material combinations.

References

- Jung R, Sailer I, Hammerle CH, et al. In-vitro color changes of soft tissues caused by restorative materials. *Int J Perio Rest Dent* 2007;27:251-257
- Park S, Da Silva J, Weber HP, et al. Optical phenomenon of peri-implant soft tissue. Part I. Spectrophotometric assessment of natural tooth gingiva and peri-implant mucosa. *Clin Oral Imp Res* 2007;18:569-574
- Manicone P, Rossi P, Raffaelli L. An overview of zirconia ceramics. Basic properties and clinical applications. *J Dent* 2007;35:819-826
- Apicella D, Veltri M, Balleri P, et al. Influence of abutment material on the fracture strength and failure modes of abutment-fixture assemblies when loaded in a bio-faithful simulation. *Clin Oral Imp Res* 2011;22:182-188
- Sailer I, Sailer T, Stawarczyk B, et al. In vitro study of the influence of the type of connection on the fracture load on zirconia abutments with internal and external implant-abutment connections. *Int J Oral Maxillofac Imp* 2009;24:850-858
- Cosgarea R, Gasparik C, Duda D, et al. Peri-implant soft tissue color around titanium and zirconia abutments. A prospective randomized controlled clinical study. *Clin Oral Imp Res* 2015;26:537-544
- Linkevicius T, Vindaxiute E, Puisys A, et al. The influence of margin location on the amount of undetected cement excess after delivery of cement-retained implant restorations. *Clin Oral Imp Res* 2011;22:1379-1384
- Korsch M, Walther W. Peri-implantitis associated with type of cement. A retrospective analysis of different types of cement and their clinical correlation to the peri-implant tissue. *Clin Imp Dent Relat Res* 2015;17:e434-e443
- Dumbrigue H, Abanomi A, Cheng L. Techniques to minimize excess luting agent in cement-retained implant restorations. *J Prosthet Dent* 2002;87:112-114
- Long L, Alqarni H, Masri R. Influence of implant abutment fabrication method on clinical outcomes: a systematic review. *Eur J Oral Implantol* 2017;10:67-77
- Wang RR, Fenton A. Titanium for prosthodontic applications. A review of the literature. *Quintessence Int* 1996;27:401-408
- Harman MK, Banks SA, Hodge WA. Wear analysis of a retrieved hip implant with titanium nitride coating. *J Arthroplast* 1997;12:938-945
- Jehn HA, Baumgärtner ME. Corrosion studies with hard coating-substrate systems. *Surface and Coat Technol* 1992;54/55:108-114
- Satomi K, Akagawa Y, Nikai H, et al. Tissue response to implanted ceramic-coated titanium alloys in rats. *J Oral Rehab* 1988;15:339-345
- Al Jabbari Y, Fehrman J, Barnes A, et al. Titanium nitride and nitrogen ion implanted coated dental materials. *Coatings* 2012;2:160-178
- Balla VK, Bhat A, Bose S, et al. Laser processed TiN reinforced Ti6Al4V composite coatings. *J Mech Behav Biomed Mater* 2012;6:9-20
- Sahasrabudhe H, Soderlind J, Bandyopadhyay A. Laser processing of *in situ* TiN/Ti composite coating on titanium. *J Mech Behav Biomed Mater* 2016;53:239-249
- Pappas MJ, Makris G, Buechel FF. Titanium nitride ceramic film against polyethylene: a 48 million cycle wear-test. *Clin Orthop* 1995;317:64-70
- Peterson CD, Hillberry BM, Heck DA. Component wear of total knee prosthesis using Ti-6Al-4V, titanium nitride coated Ti-6Al-4V and cobalt-chromium molybdenum femoral components. *J Biomed Mater Res* 1988;22:887-903
- Knotek O, Löffler F. Physical vapor deposition coatings for dental prosthesis. *Surf Coat Technol* 1992;54:536-540
- Endo K, Sachdeva R, Araki Y, et al. Effects of titanium nitride coatings on surface and corrosion characteristics of Ni-Ti alloy. *Dent Mater* 1994;13:228-239
- Mezger PR, Creugers NH. Titanium nitride coatings in clinical dentistry. *J Dent* 1992;20:342-344
- Pecnik C, Courty C, Muff D, et al. Fracture toughness of esthetic dental coating systems by nanoindentation and FIB sectional analysis. *J Mech Behav Biomed Mater* 2015;47:1-11
- Yeh TS, Wu JM, Hu LJ. The properties of TiN thin films deposited by pulsed direct current magnetron sputtering. *Thin Solid Films* 2008;516:7294-7298
- Lousa A, Esteve J, Mejia JP, et al. Influence of deposition pressure on the structural, mechanical and decorative properties of TiN thin films deposited by cathodic arc evaporation. *Vacuum* 2007;81:1507-1510
- Scarano A, Piatelli M, Vrespa G, et al. Bone Healing around titanium and titanium-nitride-coated dental implants with three surfaces: an experimental study in rats. *Clin Imp Dent Relat Res* 2003;5:103-111
- Xuereb M, Camilleri J, Attard N. Systematic review of current dental implant coating materials and novel coating techniques. *Int J Prosthodont* 2015;28:51-59
- Wang X, Bai S, Fang L, et al. Effect of plasma nitriding and titanium nitride coating on the corrosion resistance of titanium. *J Prosthet Dent* 2016;116:450-456
- Galante JO, Lemons J, Spector M, et al. The biologic effects of implant materials. *J Orthop Res* 1991;9:760-775
- Sculean A, Gruber R, Bosshardt DD. Soft tissue wound healing around teeth and dental implants. *J Clin Periodontol* 2014;41:s6-s22

31. Ricci J, Alexander H. Cell response to surfaces: a concise summary. *Int J Perio Rest Dent* 2016;36:s39-s46
32. Corvino E, Pesce P, Mura R, et al. Influence of modified titanium abutment surface on peri-implant soft-tissue behavior: a systematic review of in vitro studies. *Int J Oral Maxillofac Imp* 2020;35:503-519
33. Brunello G, Brun P, Gardin C et al. Biocompatibility and antibacterial properties of zirconium nitride coating on titanium abutments. *PLoS One* 2018;13:e0199591
34. Moon YH, Yoon MK, Moon JS, et al. Focal adhesion linker proteins expression of fibroblasts related to adhesion in response to different transmucosal abutment surfaces. *J Adv Prosthodont* 2013;5:341-350
35. Ritz U, Nusselt T, Sewing A et al. The effect of different collagen modifications for titanium and titanium nitride surfaces on functions of gingival fibroblasts. *Clin Oral Investig* 2017;21:255-265
36. Kim YS, Ko Y, Kye SB, et al. Human gingival fibroblast (HGF-1) attachment and proliferation on several abutment materials with various colors. *Int J Oral Maxillofac Imp* 2014;29:969-975
37. Kim YS, Shin SY, Moon SK, et al. Surface properties correlated with the human gingival fibroblasts attachment on various materials for implant abutments: a multiple regression analysis. *Acta Odontol Scan* 2015;73:38-47
38. Gröbner-Schreiber B, Neubert A, Muller W, et al. Fibroblast growth on surface-modified dental implants. An in-vitro study. *J Biomed Mat Res* 2003;64A:591-599
39. Gröbner-Schreiber B, Griepentrog M, Haustein I, et al. Plaque formation on surface-modified dental implants. An in-vitro study. *Clin Oral Imp Res* 2001;12:543-551
40. Gröbner-Schreiber B, Hanning M, Dück A, et al. Do different implant surfaces exposed in the oral cavity of humans show different biofilm compositions and activities?. *Europ J Oral Sci* 2004;112:516-522
41. Scarano A, Piatelli M, Vrespa G, et al. Bacterial adhesion on titanium nitride-coated and uncoated implants: an in-vivo human study. *J Oral Implantol* 2003;29:80-85
42. Quirynen M, Van der Mei HC, Bollen CM, et al. An in vivo study of the influence of the surface roughness of implants in the microbiology of supra- and subgingival plaque. *J Dent Res* 1993;72:1304-1309
43. Bollen CM, Lambrechts P, Quirynen M. Comparison of surface roughness of oral hard materials to the threshold surface roughness for bacterial plaque retention. A review of the literature. *Dent Mater* 1997;13:258-269
44. Bollen C, Papaioannou W, Van Eldere J, et al. The influence of abutment surface roughness on plaque accumulation and peri-implant mucositis. *Clin Oral Imp Res* 1996;7:201-211
45. Quirynen M, Bollen C, Papaioannou W, et al. The influence of titanium abutment surface roughness on plaque accumulation and gingivitis. Short term observations. *Int J Oral Maxillofac Imp* 1996;11:169-178
46. Park JB, Yang SM, Ko Y. Evaluation of the surface characteristics of various implant abutment materials using confocal microscopy and white light interferometry. *Imp Dent* 2015;24:650-656
47. Sanz-Sanchez I, Sanz-Martin I, Carrillo de Albornoz A, et al. Biological effect of abutment material on the stability of peri-implant marginal bone levels. A systematic review and meta-analysis. *Clin Oral Imp Res* 2018; 29(Suppl):124-144
48. Annunziata M, Oliva A, Basile MA, et al. The effect of titanium nitride-coating on the topographic and biological features of TPS implant surfaces. *J Dent* 2011;39:720-728
49. Clem WC, Konovalov VV, Chowdhury S, et al. Mesenchymal stem cell adhesion and spreading on microwave plasma-nitriding titanium alloy. *J Biomed Mater Res* 2006;76:279-287
50. Manso-Silvan M, Martinez-Duart JM, Ogueta S, et al. Development of human mesenchymal stem cells on DC sputtered titanium nitride thin films. *J Mater Sci Mater Med* 2002;13:289-293
51. Rapley JW, Swan RH, Hallmon WW, et al. The surface characteristics produced by various oral hygiene instruments and materials on titanium implant abutments. *Int J Oral Maxillofac Imp* 1990;5:47-52
52. Orton GS, Steele DL, Wolinsky LE. Dental professional's role in monitoring and maintenance of tissue-integrated prosthesis. *Int J Oral Maxillofac Imp* 1989;4:305-310
53. Tamura Y, Yokoyama A, Watari F, et al. Surface properties and biocompatibility of nitride titanium for abrasion resistant implant materials. *Dent Mater J* 2002;21:355-372
54. Sawase T, Yoshida K, Taira Y, et al. Abrasion resistance of titanium nitride coatings formed on titanium by ion-beam-assisted deposition. *J Oral Rehab* 2005;32:151-157
55. Alsabeeha N, Swain MV, Payne A. Clinical performance and material properties of single-implant overdentures attachment systems. *Int J Prosthodont* 2011;24:247-254
56. Kourtev J, Pascova R, Weibmantel E. A modified method for arc deposition of Ti-N thin films. *Vacuum* 1997;48:7-12
57. Espinoza-Beltran FJ, Che-Soberanis O, García González L, et al. Effect of the substrate bias potential on crystalline grain size, intrinsic stress and hardness of vacuum arc evaporated TiN/c-Si coatings. *Thin solid film* 2003;437:170-175
58. Stavenga DG. Thin film and multilayer optics cause structural colors of many insect and birds. *Mater Today Proceed* 2014;1S:109-121
59. Pecnik C, Roos M, Muff D, et al. In vitro color evaluation of esthetic coatings for metallic dental implants and implant prosthetic appliances. *Clin Oral Imp Res* 2015;26:563-571
60. Ishikawa-Nagai S, Da Silva JD, Weber HP, et al. Optical phenomenon of peri-implant soft tissue. Part II. Preferred implant neck color to improve soft tissue esthetics. *Clin Oral Imp Res* 2007;18:575-580
61. Mengel R, Meer C, Flores-de-Jacoby L. The treatment of uncoated and titanium nitride-coated abutments with different instruments. *Int J Oral Maxillofac Imp* 2004;19:232-238
62. Yang Y, Ong J, Tian J. Deposition of highly adhesive ZrO₂ coatings on Ti and CoCrMo implant materials using plasma nitriding. *Biomaterials* 2003;24:619-627
63. Tong Y, Gao TW, Wang J, et al. Effects of plasma nitriding and TiN coating duplex treatment on wear resistance of commercially pure titanium. *Adv Mater Res* 2011;217:1050-1055
64. Shenhar A, Gotman I, Radin S, et al. Titanium nitride-coatings on surgical titanium alloys produced by a powder immersion reaction assisted coating method: residual stresses and fretting behavior. *Surf Coat Technol* 2000;126:210-218
65. Lim HP, Lee KM, Koh YI, et al. Allergic contact stomatitis caused by a titanium nitride-coated implant abutment. A clinical report. *J Prosthet Dent* 2012;108:209-213
66. Ferrari M, Tricarico MG, Cagidiaco MC, et al. 3-year randomized controlled prospective clinical trial on different CAD-CAM implant abutments. *Clin Imp Dent Relat Res* 2016;18:1134-1141
67. Parpaiola A, Toia M, Norton M, et al. CAD/CAM implant abutments: peri-implant hard and soft tissue response with up to 4 years of follow-up. A retrospective cohort study evaluation. *Int J Period Rest Dent* 2020;40:193-201

68. Ferrari M, Cagidiaco MC, García-Godoy K, et al. Effect of different prosthetic abutments on peri-implant soft tissue. A randomized controlled clinical trial. *Am J Dent* 2015;28:85-89
69. Borges T, Lima T, Carvalho A, et al. Clinical outcome of inter-proximal papilla between a tooth and a single implant treated with CAD/CAM abutments: a cross-sectional study. *J Oral Maxillofac Res* 2012;3:e4
70. Borges T, Lima T, Carvalho A, et al. The influence of customized abutments and custom metal abutments on the presence of the interproximal papilla at implants inserted in single-unit gaps: a 1-year prospective clinical trial. *Clin Oral Imp Res* 2014;25:1222-1227
71. Happe A, Schulte-Mattler V, Fickl S, et al. Spectrophotometric assessment of peri-implant mucosa after restoration with zirconia abutments veneered with fluorescent ceramic: a controlled, retrospective clinical study. *Clin Oral Imp Res* 2013;24:28-33
72. Van Brakel R, Noordmans HJ, Frenken J, et al. The effect of zirconia and titanium implant abutments on light reflection of the supporting soft tissues. *Clin Oral Imp Res* 2011;22:1173-1178
73. Bressan E, Paniz G, Lops D, et al. Influence of abutment material on gingival color of implant-supported all-ceramic restorations: a prospective multicenter study. *Clin Oral Imp Res* 2011;22:631-637
74. Kim A, Campbell SD, Viana MA, et al. Abutment material effect on peri-implant soft-tissue color and perceived esthetics. *J Prosthodont* 2016;25:634-640
75. Lops D, Stellini E, Sbricoli L, et al. Influence of abutment material on peri-implant soft tissues in anterior areas with thin gingival biotype. A multicentric prospective study. *Clin Oral Imp Res* 2017;28:1263-1268
76. Abduo J, Lyons K. Rationale for the use of CAD/CAM technology in implant prosthodontics. *Int J Dent* 2013;2013:768121
77. Wasiluk G, Chomik E, Gehrke P, et al. Incidence of undetected cement on CAD-CAM monolithic zirconia crown and customized CAD-CAM implant abutments. A prospective case series. *Clin Oral Imp Res* 2017;28:774-778
78. Ramsey CD, Leal KF. Patient-specific abutments for anterior maxillary implants. *Comp Contin Educ Dent* 2014;35:404-411
79. Fabbri G, Fraedani M, Dellificorelli G, et al. Clinical evaluation of the influence of connection type and restoration height on the reliability of zirconia abutments: a retrospective study on 965 abutments with mean 6-year follow-up. *Int J Perio Rest Dent* 2017;37:19-31
80. Sailer I, Asgeirsson A, Thoma D, et al. Fracture strength of zirconia implant abutments on narrow diameter implants with internal and external implant abutment connections. A study on the titanium resin base concept. *Clin Oral Imp Res* 2018;29:411-423
81. Massoud SN, Hunter J, Holdsworth B, et al. Early femoral loosening in one design of cemented hip replacement. *J Bone Joint Surg* 1997;79-B:603-608
82. Lombardi AV, Mallory TH, Vaughan BK, et al. Aseptic loosening in total hip arthroplasty secondary to osteolysis induced by wear debris from titanium-alloy modular femoral heads. *J Bone Joint Surg* 1989;71-A:1337-1342
83. Shamseer L, Moher D, Clarke M, et al. Preferred reporting items for systematic review and meta-analysis protocols (PRISMA-P) 2015: elaboration and explanation. *BMJ* 2015;350:g7647